
pekatherm

warm softness

- ES** INSTRUCCIONES DE USO
VENDAJES TÉRMICOS
- EN** OPERATING INSTRUCTIONS
HEATING BANDAGES
- RU** ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
БАНДАЖИ С ПОДОГРЕВОМ
- PT** INSTRUÇÕES DE UTILIZAÇÃO
BANDAGENS AQUECEDORAS
- FR** MANUEL D'UTILISATION
BANDAGES CHAUFFANTS SANS FIL
- DE** GEBRAUCHSANLEITUNG
THERMOBANDAGEN
- IT** ISTRUZIONI D'USO
SEGMENTI RISCALDANTI
- HE** הוראות שימוש
מגינים חימום

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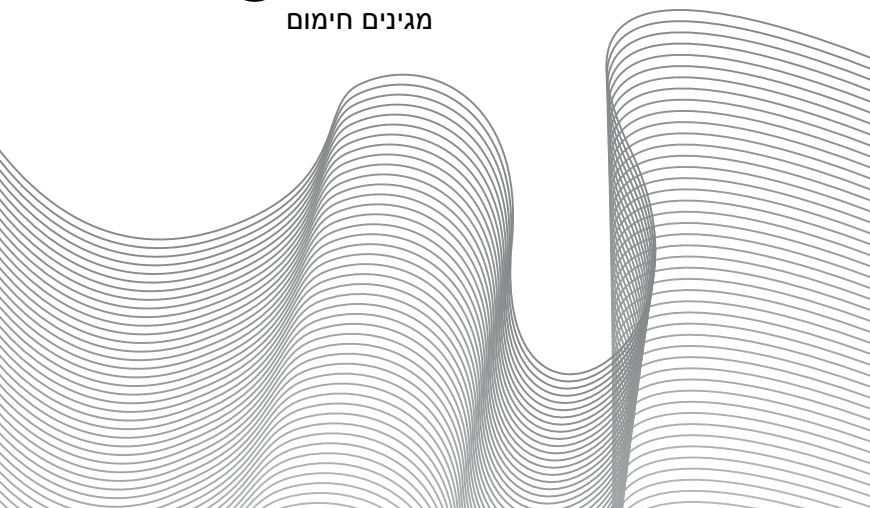IMPORTANT INFORMATION! PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE!

THESE INSTRUCTIONS PROVIDE INFORMATION AND ADVICE ON INSTALLATION, USE AND MAINTENANCE OF THE PRODUCT. READ THEM CAREFULLY BEFORE USING THE PRODUCT TO ENSURE SAFE AND PROPER USE. NON OBSERVANCE OF THESE INSTRUCTIONS CAN CAUSE PERSONAL INJURIES, AS WELL AS PRODUCT DAMAGE OR DESTRUCTION. KEEP THESE INSTRUCTIONS AND PASS THEM ON WITH THE PRODUCT IF SOLD OR LENT TO OTHERS.

PRODUCT PARTS

(see pic. 1, 2, 3 and 4)

This set comprises heating bandage, battery, battery pocket, charger and user manual.

- | | |
|---|--------------------------|
| a. Heating bandage (according to model) | f. Bandage connector |
| b. Battery | g. Battery connector |
| c. Battery pocket | h. Charger LED indicator |
| d. Charger | i. Charger connector |
| e. User manual | |

BEFORE USE

1. After removing the packaging, inspect the product carefully to make sure it has no signs of damage caused during transport. In case of any doubt, do not use the product and take it to your dealer or to an authorized service centre.
2. Use only the battery and the charger supplied with the product. Never interchange them with those of other appliances.
3. Make sure your mains rating corresponds to the stated on the charger's label.
4. This product must not be used with sick or disabled persons, with persons insensitive to heat or persons with restricted physical, sensory or mental capabilities. It also must not be used by children or by people lacking the required experience and knowledge of product operation unless they are properly instructed on how to operate the product safely. In all circumstances, while using the product, children should always be supervised by an adult to prevent them from playing with it.
5. Do not use the product under the effects of alcohol, sleeping pills or sedatives.
6. Do not use the product on animals.
7. Do not apply the product to the body areas where the skin presents scars, wounds, inflammation or any other alteration due to illness. The product should not be used by persons suffering from phlebitis, arthritis or other inflammations since heat may aggravate the disease. Consult your physician.
8. This product is designed for home use only. It is not intended for professional applications or medical use in hospitals.
9. Heating bandages should only be used for warming a specific part of human body. Any other use of the product will be considered hazardous. The manufacturer will decline all responsibility for any damage due to improper or incorrect use.
10. Pacemaker users should consult their physician before using the product.
11. Do not use the product folded or creased. Do not store the product creased.
12. Do not stick pins, needles or other sharp objects into the product.
13. Do not place heavy objects onto the product.
14. Never bring the product into contact with other sources of heat, such as other heating bandages, heating pads, heating underblankets, hot water bottles, etc.

15. Do not cover the battery while operating the product; do not place it on the heating element or under it while the unit is in operation.
16. Excessively long use of heating bandages running at a high temperature may produce skin burns. Avoid falling asleep with heating bandage in operation.
17. It is not recommended to apply the bandage directly on the skin for persons with sensitive skin. Take care to inspect frequently the skin in order to avoid skin burns. Disconnect the bandage immediately if you note burning or if the skin turns red or presents inflammations.
18. Do not leave the product unattended while in operation.
19. In any circumstances avoid wetting the battery and the charger. Do not use the product if wet.
20. Do not expose the product to environmental agents (rain, sun, etc.).
21. Keep the bandage, the battery and the charger away from fire.
22. Do not use creams or body oils that could damage the bandage.
23. Detach the battery and unplug the charger while not in use.
24. Detach the battery before cleaning the bandage.
25. For bandage cleaning please follow instructions in a specific section of this manual.
26. Do not hold the product by the battery to move or carry it.
27. Inspect the product frequently for signs of wear or damage. If these signs appear, do not use the product and take it to an authorized service centre.
28. This product has no parts intended to be repaired or replaced by the user. If any part of the product fails, or the product does not perform properly, or if there is any other kind of anomaly, do not attempt to repair it and instead take the product to an authorized service centre.
29. For any further information on our products please consult our customer service.

ADVICE ON BATTERY USE

30. Do not use the battery with other appliances.
31. Do not use the battery if the connector is damaged or presents signs of rust.
32. Never short-circuit the battery. Avoid keeping it together with small metallic objects. Do not wrap it in any conductor material.
33. Do not attempt to disassemble or alter any part of the battery.
34. Keep the battery away from children.
35. Keep the battery away from fire and water.
36. Battery is supplied half charged. Proceed to fully charge it before use.
37. To recharge the battery use only the supplied charger.
38. Once the battery is fully charged unplug it from the charger. Do not leave it plugged in when unnecessary.
39. Batteries gradually self-discharge with time. Proceed to fully charge the battery before each use.
40. Before storing the battery for a long period of time half charge it as explained in a specific section of this manual. Repeat the operation at least every half a year.
41. The capacity of rechargeable batteries gradually decreases. The supplied battery is designed to provide 300 discharge /charge cycles.
42. After a battery has reached the end of its lifetime take it to a green collecting point in order to assure the subsequent dispatch of the discarded elements for recycling, treatment and environmentally compatible disposal. Do not dispose it in fire or in water.

HEATING BANDAGES

Our heating bandages are developed and created in collaboration with physiotherapists and constitute the ultimate alternative in thermal therapy:

- permit local targeted heat treatment. All our bandages are uniquely designed to follow the complex shape of a specific body region.
- neoprene-based multilayer fabric structure, highly elastic material permits perfect fitting and light compression on the joint where the bandage is applied.

INSTRUCTIONS

Insulating material retains heat and maintains the joint and adjacent muscles hot during treatment. Thermal insulation helps to prevent heat dissipation concentrating the heat in the affected area, which contributes to an easier heat penetration into the joint.

- low voltage operation assures maximum safety of use.
- its battery powered technology provides unlimited motion to the user:

These products will actually make you forget the well known traditional heating products like hot water bottles, gel packs needing microwave or boiling water activating, etc.

- Technical characteristics:

Material	Neoprene
Heating element	Resistive wire
Power supply	Battery
Battery characteristics	Li-ion 3.7V 2200mAh
Power	12W
Max. temperature, aprox.	55°C [‡]
Battery autonomy	150min [‡]
Charger	Input 100-240VAC 50-60Hz, 4.2V 1.0A
Charging time, aprox.	4h

[‡] Note: all measures were collected in specific laboratory conditions

BENEFITS OF HEAT THERAPY

Local heat application stimulates metabolism and cell functions, dilates blood vessels increasing blood flow, oxygenation and nutrition of the affected area. This contributes to stimulate processes of regeneration and recovery of affected tissues, as well as elimination of toxins, products of tissue decay and lactic acid which improves vitality and helps to prevent other diseases.

Heat also reduces muscular tension and stiffness caused by overuse or overstretching and is also very helpful in case of rheumatic pains. Heat also relieves pain sensation since it reduces pain receptors stimulation and transmitting pain signals to the spinal cord and brain. Heat eases muscle tension and stress, and definitely makes great improvement in our sensation of comfort and well-being.

REASONABLE HEAT TREATMENT

In order to bring positive effects, heat treatment does not require high temperatures or long time applications. It is with this premise in mind that the new series of heating bandages was created. Operating temperature is set at approximately 55°C, as measured in laboratory conditions. To get maximum benefit from heat treatment it is recommended to apply the bandage during 20 to 30 minutes. Consult your physician or physiotherapist.

Heat must be clearly perceptible but never excessive. It should always be pleasant and should never cause physical trouble. Take care to inspect frequently the skin in order to avoid skin burns. Disconnect the bandage immediately if you note burning or itch or if the skin turns red or presents inflammations.

USE INSTRUCTIONS

1. Put the bandage onto a specific body part and fix it with the help of elastic straps.
2. Connect the battery to the bandage, and check that the battery has been properly connected.
3. Insert the battery in the battery pocket, (fig. 4) and stick it on any area of the bandage which is most comfortable.
4. At this point the bandage begins warming up and in a few seconds you will be able to perceive a pleasant heat. Make sure the bandage is fixed correctly. Check that there are no creases or folds.
5. You may extend the treatment period if necessary; remember, however, that generally it should not last more than 20 to 30 minutes (see "Reasonable heat treatment").
6. If the battery is not to be used for a considerable period of time store it disconnected from the bandage. Half charge the battery before storing as explained in a specific section of this manual.

BATTERY CHARGING

1. Connect the battery to the charger and plug the charger into the mains socket. (fig. 3) Red LED indicates the charging cycle is on.

NOTE: Green LED indicates that the battery is already fully charged. In this case unplug it from the charger and unplug the charger from the mains socket.

2. Charging cycle will take about 4 hours for a completely discharged battery. When charging is over LED turns green. When the battery is fully charged disconnect it from the charger. Unplug the charger from the mains to reduce consumption.

NOTE: Half charge the battery before storing if it is not to be used for a considerable period of time. Proceed as follows: completely discharge the battery connecting it to the bandage and then put it to charge for about 2 hours. Repeat the procedure at least once in half a year.

TROUBLESHOOTING TABLE

PROBLEM	SOLUTION
1. The bandage does not warm up	The battery may be connected to the bandage incorrectly. Check the connection. Remember that the battery may be completely discharged. Disconnect the battery from the bandage, and proceed to fully charge (see chapter on battery charging). It may also be a question of feeling. Remember that quick checks carried out with uncovered bandage are not reliable. Try to cover the product with a blanket or a towel and wait at least 30 minutes. If you still do not perceive heat, take the product to an authorized service centre.
2. The bandage warms up insufficiently	If you perceive that the heat emitted by the product is insufficient, try the following: remove the clothes from the part of the body to which the bandage is intended to be applied and then apply the bandage directly on the skin, covering its other side with a towel. This will increase heat transmission, as well as avoid the dissipation of some of the heat.
3. Battery does not charge.	Check the connection between the battery and the charger, and check if the charger is plugged in correctly. Check if your mains socket has voltage. If the problem persists take the product to an authorized service centre.

INSTRUCTIONS

CLEANING

Before cleaning the bandage disconnect the battery. Never wet the battery and the charger. For a basic cleaning of the bandages soft brushing is recommended. To get smaller stains out use a moist cloth or sponge with some liquid soap or cleaning foam.

Apply the foam regularly on the product's surface and then brush it. If washing is necessary remember that only hand washing is permitted. Use a soft detergent and rub or brush softly avoiding wringing the bandage. Never bleach the bandage.

After washing, extend the bandage flat or hang it, avoid folding it. Always place the bandage in the shadow since direct sunlight could damage the bandage. Do not iron the bandage.

STORAGE

Let the product cool down before folding it for storage. Avoid excessive folds or creases. Store it in a cool dry place. Do not place any heavy objects on the product. Follow the instructions on battery storage when not in use for long periods of time.

Inspect the product (the bandage, the battery and the charger) frequently for signs of wear, damage or misuse. If these signs appear, do not use the product. Remember this product has no parts intended to be repaired or replaced by the user. Do not try to carry out repairs and, instead, take the product to an authorized service centre.

ADVICE ON DISPOSAL

According to EU Directive 2012/19/EC WEEE (Waste of Electrical and Electronic Equipment), electrical and electronic devices must be disposed of separately from household refuse, delivering them to a special refuse collection point, in order to assure the subsequent dispatch of the discarded equipment for recycling, treatment and environmentally compatible disposal, thus contributing to the prevention of possible negative effects on the environment and on health. Your local authorities will provide you further information on how to dispose of your inoperative appliances.

WARRANTY

pekatherm S.L. guarantees this product against manufacturing and material defects for a period of two years from the date of purchase.

pekatherm S.L. will replace or repair the defective item free of charge, including labour and spare parts, if the defect derives from manufacturing and/or material reasons, and if the product is handed by the customer to an authorized service centre.

If replacement or repair is not possible, the customer may claim for refund or buy a superior model paying the difference. This warranty will not be valid if not accompanied by the original purchase bill or ticket. Application of the warranty will not imply any extension of its period. This warranty will not cover any defect due to reasons beyond manufacturing, such as transportation damages, non observance of the operating instructions, careless or improper use, wear and tear, damages caused by a blow or fall, water, chemical agents, etc.

This warranty will not cover any cases of improper use of the product, such as professional or medical use. This warranty will not cover any damage if improper repairs were attempted by unauthorized personnel.

pekatherm S.L. declines all responsibility for any damage caused directly or indirectly to persons, pets or properties deriving from non observance of these operating instructions. This warranty is only valid in the country where the product was purchased.

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